

SCD1(M38) Antibody

✓ 100 µl
(10 Western mini-blot)

Orders ■ 877-616-CELL (2355)
orders@cellsignal.com

Support ■ 877-678-TECH (8324)
info@cellsignal.com

Web ■ www.cellsignal.com

rev. 10/22/08

This product is for *in vitro* research use only and is not intended for use in humans or animals.
This product is not intended for use as a therapeutic or in diagnostic procedures.


Entrez-Gene ID #6319
Swiss-Prot Acc. #O00767

Applications	Species Cross-Reactivity*	Molecular Wt.	Source
W, IF-IC Endogenous	M, (H)	37 kDa	Rabbit**

Background: Stearoyl-CoA desaturase 1 (SCD1), a key lipogenic enzyme in the endoplasmic reticulum, catalyzes the conversion of palmitoyl-CoA and stearoyl-CoA to palmitoleoyl-CoA (16:1) and oleoyl-CoA (18:1) (1,2,3). Palmitoleate and oleate are the major components of triglycerides, membrane phospholipids and cholesterol esters (1). The SCD1-knockout mice show improved insulin sensitivity and reduced body fat (1). Disruption of SCD1 in mouse brown adipose tissue strengthens insulin signaling, resulting in the increased translocation of Glut4 to plasma membrane and therefore enhanced uptake of glucose into brown adipocytes (4). Furthermore, SCD1 was shown to be essential for the onset of diet-induced body weight gain (1) and insulin resistance in liver (5).

Specificity/Sensitivity: SCD1 (M38) Antibody detects endogenous levels of total SCD1 protein.

Source/Purification: Polyclonal antibodies are produced by immunizing animals with a synthetic peptide (KLH-coupled) derived from the sequence around M38 of human SCD1. Antibodies are purified by protein A and peptide affinity chromatography.

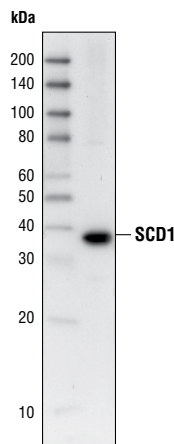

Western blot analysis of 3T3-L1 cell lysates, using SCD1 (M38) Antibody.

Background References:

- (1) Ntambi, J.M. et al. (2002) *Proc. Natl. Acad. Sci. USA* 99, 11482–11486.
- (2) Kato, H. et al. (2006) *J. Cell. Sci.* 119, 2342–2353.
- (3) Ozols, J. (1997) *Mol. Biol. Cell* 8, 2281–2290.
- (4) Rahman, S.M. et al. (2005) *Am. J. Physiol. Endocrinol. Metab.* 288, E381–387.
- (5) Gutiérrez-Juárez, R. et al. (2006) *J. Clin. Invest.* 116, 1686–1695.

Storage: Supplied in 10 mM sodium HEPES (pH 7.5), 150 mM NaCl, 100 µg/ml BSA and 50% glycerol. Store at -20°C.

Do not aliquot the antibody.

*Species cross-reactivity is determined by Western blot.


**Anti-rabbit secondary antibodies must be used to detect this antibody.

Recommended Antibody Dilutions:

Western blotting 1:1000
Immunofluorescence (IF-IC) 1:50

Companion Products:

SCD1 (R347) Antibody #2283
Fatty Acid Synthase Antibody #3189
ATGL Antibody #2138
Phospho-HSL (Ser563) Antibody #4139
Phospho-HSL (Ser565) Antibody #4137
Phospho-HSL (Ser660) Antibody #4126
Phototope®-HRP Western Blot Detection System, Anti-rabbit IgG, HRP-linked Antibody #7071
Anti-rabbit IgG, HRP-linked Antibody #7074
Prestained Protein Marker, Broad Range (Premixed Format) #7720
Biotinylated Protein Ladder Detection Pack #7727
20X LumiGLO® Reagent and 20X Peroxide #7003

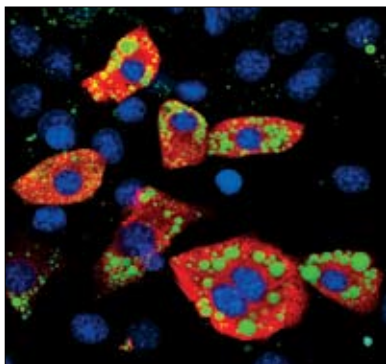

Confocal immunofluorescent analysis of 3T3-L1 cells, using SCD1 (M38) Antibody (red) showing cytoplasmic localization in differentiated cells. Lipid droplets have been labeled with BODIPY 493/503 (green). Blue pseudocolor = DRAQ5™ (fluorescent DNA dye).

IMPORTANT: For Western blots, incubate membrane with diluted antibody in 5% w/v BSA, 1X TBS, 0.1% Tween-20 at 4°C with gentle shaking, overnight.

Applications Key: W—Western IP—Immunoprecipitation IHC—Immunohistochemistry ChIP—Chromatin Immunoprecipitation IF—Immunofluorescence F—Flow cytometry E-P—ELISA-Peptide
Species Cross-Reactivity Key: H—human M—mouse R—rat Hm—hamster Mk—monkey Mi—mink C—chicken Dm—D. melanogaster X—Xenopus Z—zebra fish B—bovine
Dg—dog Pg—pig Sc—S. cerevisiae All—all species expected Species enclosed in parentheses are predicted to react based on 100% sequence homology.

Western Immunoblotting Protocol (Primary Antibody Incubation in BSA)

For Western blots, incubate membrane with diluted antibody in 5% w/v BSA, 1X TBS, 0.1% Tween-20 at 4°C with gentle shaking, overnight.

A Solutions and Reagents

NOTE: Prepare solutions with Milli-Q or equivalently purified water.

- 1X Phosphate Buffered Saline (PBS)
- 1X SDS Sample Buffer:** 62.5 mM Tris-HCl (pH 6.8 at 25°C), 2% w/v SDS, 10% glycerol, 50 mM DTT, 0.01% w/v bromophenol blue or phenol red
- Transfer Buffer:** 25 mM Tris base, 0.2 M glycine, 20% methanol (pH 8.5)
- 10X Tris Buffered Saline (TBS):** To prepare 1 liter of 10X TBS: 24.2 g Tris base, 80 g NaCl; adjust pH to 7.6 with HCl (use at 1X).
- Nonfat Dry Milk (weight to volume [w/v])
- Blocking Buffer:** 1X TBS, 0.1% Tween-20 with 5% w/v nonfat dry milk; for 150 ml, add 15 ml 10X TBS to 135 ml water, mix. Add 7.5 g nonfat dry milk and mix well. While stirring, add 0.15 ml Tween-20 (100%).
- Wash Buffer:** 1X TBS, 0.1% Tween-20 (TBS/T)
- Bovine Serum Albumin (BSA)
- Primary Antibody Dilution Buffer:** 1X TBS, 0.1% Tween-20 with 5% BSA; for 20 ml, add 2 ml 10X TBS to 18 ml water, mix. Add 1.0 g BSA and mix well. While stirring, add 20 µl Tween-20 (100%).
- Phototope[®]-HRP Western Blot Detection System #7071:** Includes biotinylated protein ladder, secondary anti-rabbit (#7074) antibody conjugated to horseradish peroxidase (HRP), anti-biotin antibody conjugated to HRP, LumiGLO[®] chemiluminescent reagent and peroxide.
- Prestained Protein Marker, Broad Range (Premixed Format) #7720
- Biotinylated Protein Ladder Detection Pack #7727
- Blotting Membrane:** This protocol has been optimized for nitrocellulose membranes, which CST recommends. PVDF membranes may also be used.

B Protein Blotting

A general protocol for sample preparation is described below.

- Treat cells by adding fresh media containing regulator for desired time.
- Aspirate media from cultures; wash cells with 1X PBS; aspirate.
- Lyse cells by adding 1X SDS sample buffer (100 µl per well of 6-well plate or 500 µl per plate of 10 cm diameter plate). Immediately scrape the cells off the plate and transfer the extract to a microcentrifuge tube. Keep on ice.
- Sonicate for 10–15 seconds to shear DNA and reduce sample viscosity.
- Heat a 20 µl sample to 95–100°C for 5 minutes; cool on ice.
- Microcentrifuge for 5 minutes.
- Load 20 µl onto SDS-PAGE gel (10 cm x 10 cm).

NOTE: CST recommends loading prestained molecular weight markers (#7720, 10 µl/lane) to verify electrotransfer and biotinylated protein ladder (#7727, 10 µl/lane) to determine molecular weights.

- Electrotransfer to nitrocellulose or PVDF membrane.

C Membrane Blocking and Antibody Incubations

NOTE: Volumes are for 10 cm x 10 cm (100 cm²) of membrane; for different sized membranes, adjust volumes accordingly.

- (Optional) After transfer, wash nitrocellulose membrane with 25 ml TBS for 5 minutes at room temperature.
- Incubate membrane in 25 ml of blocking buffer for 1 hour at room temperature.
- Wash three times for 5 minutes each with 15 ml of TBS/T.
- Incubate membrane and primary antibody (at the appropriate dilution) in 10 ml primary antibody dilution buffer with gentle agitation overnight at 4°C.
- Wash three times for 5 minutes each with 15 ml of TBS/T.
- Incubate membrane with HRP-conjugated secondary antibody (1:2000) and HRP-conjugated anti-biotin antibody (1:1000) to detect biotinylated protein markers in 10 ml of blocking buffer with gentle agitation for 1 hour at room temperature.
- Wash three times for 5 minutes each with 15 ml of TBS/T.

D Detection of Proteins

- Incubate membrane with 10 ml LumiGLO[®] (0.5 ml 20X LumiGLO[®], 0.5 ml 20X Peroxide and 9.0 ml Milli-Q water) with gentle agitation for 1 minute at room temperature.

NOTE: LumiGLO[®] substrate can be further diluted if signal response is too fast.

- Drain membrane of excess developing solution (do not let dry), wrap in plastic wrap and expose to x-ray film. An initial 10-second exposure should indicate the proper exposure time.

NOTE: Due to the kinetics of the detection reaction, signal is most intense immediately following LumiGLO[®] incubation and declines over the following 2 hours.

Immunofluorescence Protocol

***IMPORTANT:** Please refer to the **APPLICATIONS** section on the front page of the data sheet to determine **IF THIS PRODUCT** is validated and approved for the specific protocol you will be using.

A Solutions and Reagents

NOTE: Prepare solutions with Milli-Q or equivalently purified water.

- 1. 10X Phosphate Buffered Saline (PBS):** To prepare 1 L add 80 g sodium chloride (NaCl), 2 g potassium chloride (KCl), 14.4 g sodium phosphate, dibasic (Na_2HPO_4) and 2.4 g potassium phosphate, monobasic (KH_2PO_4) to 1 L dH_2O . Adjust pH to 7.4.
- Formaldehyde, 16%, methanol free, Polysciences, Inc. (cat# 18814), use fresh, store opened vials at 4°C in dark, dilute in PBS for use.
- Xylene
- Ethanol, anhydrous denatured, histological grade, 100% and 95%
- Distilled water (dH_2O)
- Blocking Buffer:** To prepare 25 mL, add 2.5 mL 10X PBS, 1.25 mL normal serum from the same species as the secondary antibody (eg. normal goat serum, normal donkey serum) and 21.25 mL dH_2O and mix well. While stirring, add 75 μL Triton X-100 (100%).
- Antibody Dilution Buffer:** To prepare 40 mL, add 4 mL 10X PBS to 36 mL dH_2O , mix. Add 0.4 g BSA and mix well. While stirring, add 120 μL Triton X-100 (100%).
- 10 mM Sodium Citrate Buffer:** To prepare 1 L, add 2.94 g sodium citrate trisodium salt dihydrate ($\text{C}_6\text{H}_5\text{Na}_3\text{O}_7 \cdot 2\text{H}_2\text{O}$) to 1 L dH_2O . Adjust pH to 6.0.
- 1X PBS, high salt (0.4M) (high salt PBS):** To prepare 1L, add 100 ml 10X PBS to 900 ml dH_2O . Add 23.38 g NaCl and mix.
- Fluorochrome-conjugated secondary antibody

NOTE: When using any primary or fluorochrome-conjugated secondary antibody for the first time, titrate the antibody to determine which dilution allows for the strongest specific signal with the least background for your sample.

- Prolong® Gold Antifade Reagent (Invitrogen, Eugene, OR, Cat# P36930)

B Specimen Preparation

I. Cultured Cell Lines (IF-IC)

IMPORTANT: Please check the **APPLICATIONS** section of the data sheet to verify that this product is validated and approved for **(IF-IC)**.

NOTE: Cells should be grown, treated, fixed, and stained directly in multiwell plates, chamber slides, or on coverslips.

- Rinse cells briefly in PBS.
- Aspirate PBS, cover cells to a depth of 2-3 mm with 2-4% formaldehyde in PBS.

NOTE: Formaldehyde is toxic, use only in fume hood.

- Allow cells to fix for 15 minutes at room temperature.
- Aspirate fixative, rinse three times in PBS for 5 minutes each.
- Methanol Permeabilization Step (if required, please refer to front page):** After formaldehyde fixation, cover cells with ice-cold 100% methanol (use enough to cover cells completely to a depth of 3-5 mm, DO NOT LET CELLS DRY), incubate cells in methanol for 10 minutes at -20°C , rinse in PBS for 5 minutes.
- Proceed with Immunostaining section C.

II. Paraffin Sections (IF-P)

IMPORTANT: Please check the **APPLICATIONS** section of the data sheet to verify that this product is validated and approved for **(IF-P)**.

Deparaffinization/Rehydration:

- Incubate sections in three washes of xylene for 5 minutes each.
- Incubate sections in two washes of 100% ethanol for 10 minutes each.
- Incubate sections in two washes of 95% ethanol for 10 minutes each.
- Rinse sections twice in dH_2O for 5 minutes each.

Antigen Unmasking:

- Place slides in room temperature 10 mM sodium citrate buffer pH 6.0.
- Bring slides to boiling in sodium citrate buffer using water bath or microwave, then maintain at $95-99^\circ\text{C}$ for 10 minutes.
- Cool slides for 30 minutes on bench top.
- Rinse sections in dH_2O three times for 5 minutes each.
- Rinse sections in PBS for 5 minutes.
- Proceed with Immunostaining section C.

III. Frozen/Cryostat Sections (IF-F)

IMPORTANT: Please check the **APPLICATIONS** section of the data sheet to verify that this product is validated and approved for **(IF-F)**.

NOTE: Fresh frozen/unfixed sections should be fixed immediately in 2-4% formaldehyde as follows to preserve signaling epitopes.

- Cover sections with 2-4% formaldehyde in PBS

NOTE: Formaldehyde is toxic, use only in fume hood.

- Allow sections to fix for 15 minutes at room temperature.
- Rinse slides three times in PBS for 5 minutes each.

C Immunostaining

NOTE: All subsequent incubations should be carried out at room temperature unless otherwise noted in a humid light-tight box or covered dish/plate to prevent drying and fluorochrome fading.

- Block specimen in Blocking Buffer for 60 minutes.
- While blocking, prepare primary antibody by diluting as indicated on datasheet in Antibody Dilution Buffer.
- Aspirate blocking solution, apply diluted primary antibody.

NOTE: For double-labeling, prepare a cocktail of the primary antibodies at their appropriate dilution in Antibody Dilution Buffer.

- Incubate overnight at 4°C .
- Rinse three times in PBS for 5 minutes each.

OPTION: To decrease background stain, rinse in high salt PBS for two minutes between second and third PBS rinses. Be aware, this may reduce specific staining of some antibodies.

NOTE: If using primary antibodies directly conjugated with Alexa Fluor® fluorochromes, then skip to step C8.

- Incubate specimen in fluorochrome-conjugated secondary antibody diluted in Antibody Dilution Buffer for 1-2 hours at room temperature in dark.

NOTE: For double-labeling, prepare a cocktail of fluorochrome-conjugated secondary antibodies at their appropriate dilutions in Antibody Dilution Buffer.

- Rinse in PBS/high salt PBS as in step 5.
- Coverslip slides with Prolong® Gold Antifade Reagent or apply just enough to cover cells in multiwell plate.
- Seal slides by painting around edges of coverslips with nail polish.
- For best results examine specimens immediately using appropriate excitation wavelength. For long term storage, store slides flat at 4°C protected from light.