

#2761 Store at -20°C

Numb (C44B4) Rabbit mAb

✓ 100 µl
(10 Western mini-blot)



Orders ■ 877-616-CELL (2355)
orders@cellsignal.com
Support ■ 877-678-TECH (8324)
info@cellsignal.com
Web ■ www.cellsignal.com

New 10/07


This product is for *in vitro* research use only and is not intended for use in humans or animals.

Applications	Species Cross-Reactivity*	Molecular Wt.	Source	Isotype
W, IP, IF-IC, F	H	72, 74 kDa	Rabbit**	IgG


Background: Numb contains an amino-terminal phosphotyrosine-binding (PTB) domain and carboxy-terminal endocytic binding motifs for α -adaptin and EH (Eps15 homology) domain-containing proteins, indicating a role in endocytosis (1,2). There are four mammalian Numb splicing isoforms that are differentially expressed and may have distinct functions (3-5). Numb acts as a negative regulator of Notch signaling by promoting ubiquitination and degradation of Notch (6). The protein is asymmetrically segregated into one daughter cell during cell division; therefore the two daughter cells have different responses to Notch signaling and have different cell fates (7,8). The localization of Numb can also be regulated by G-protein coupled receptor (GPCR) and PKC signaling (9).

Specificity/Sensitivity: Numb (C44B4) Rabbit mAb detects endogenous levels of Numb protein.


Source/Purification: Monoclonal antibody is generated by immunizing rabbits with a synthetic peptide (Blue Carrier-coupled) corresponding to residues surrounding Phe482 of human Numb protein.

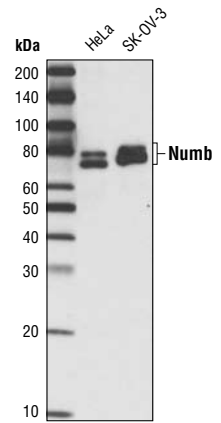

Western blot analysis of extracts from HeLa and SK-OV-3 cells using Numb (C44B4) Rabbit mAb.

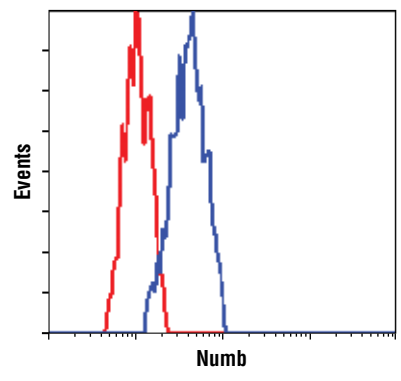

Flow cytometric analysis of A204 cells using Numb (C44B4) Rabbit mAb (blue) compared to a nonspecific negative control antibody (red).

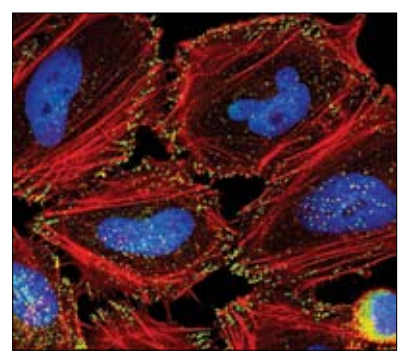

Confocal immunofluorescent analysis of HeLa cells using Numb (C44B4) Rabbit mAb (green). Actin filaments have been labeled with DY-554 phalloidin (red). Blue pseudocolor = DRAQ5™ (fluorescent DNA dye).

Entrez-Gene ID #8650
Swiss-Prot Acc. #P49757

Storage: Supplied in 10 mM sodium HEPES (pH 7.5), 150 mM NaCl, 100 µg/ml BSA, 50% glycerol and less than 0.02% sodium azide. Store at -20°C. Do not aliquot the antibody.

*Species cross-reactivity is determined by Western blot.

**Anti-rabbit secondary antibodies must be used to detect this antibody.

Recommended Antibody Dilutions:

Western Blotting	1:1000
Immunoprecipitation	1:50
Immunofluorescence (IF-IC)	1:100
Flow Cytometry	1:200

Companion Products:

- Cleaved Notch 1 (Val1744) Antibody #2421
- Notch1 (6A5) Rabbit mAb #2495
- Notch2 (8A1) Rabbit mAb #2420
- Notch4 (L5C5) Mouse mAb #2423
- DLL1 Antibody #2588
- DLL3 (G93) Antibody #2483
- DLL4 Antibody #2589
- Jagged1 (1C4) Rabbit mAb #2155
- Jagged1 (28H8) Rabbit mAb #2620
- Jagged2 (C23D2) Rabbit mAb #2210
- Jagged2 (C83A8) Rabbit mAb #2205
- MAML1 Antibody #4608
- MAML2 Antibody #4618
- Phototope®-HRP Western Blot Detection System, Anti-rabbit IgG, HRP-linked Antibody #7071
- Anti-rabbit IgG, HRP-linked Antibody #7074
- Prestained Protein Marker, Broad Range (Premixed Format) #7720

Background References:

- (1) Berdnik, D. et al. (2002) *Dev. Cell* 3, 221-231.
- (2) Santolini, E. et al. (2000) *J. Cell Biol.* 151, 1345-1352.
- (3) Dho, S.E. et al. (1999) *J. Biol. Chem.* 274, 33097-33104.
- (4) Verdi, J.M. et al. (1999) *Proc. Natl. Acad. Sci. USA* 96, 10472-10476.
- (5) Toriya, M. et al. (2006) *Dev. Neurosci.* 28, 142-155.
- (6) McGill, M.A. and McGlade, C.J. (2003) *J. Biol. Chem.* 278, 23196-23203.
- (7) Verdi, J.M. et al. (1996) *Curr. Biol.* 6, 1134-1145.
- (8) Reugels, A.M. et al. (2006) *Dev. Dyn.* 235, 934-948.
- (9) Dho, S.E. et al. (2006) *Mol. Biol. Cell.*, Epub Ahead of Print.

IMPORTANT: For Western blots, incubate membrane with diluted antibody in 5% w/v BSA, 1X TBS, 0.1% Tween-20 at 4°C with gentle shaking, overnight.

Applications Key: W—Western IP—Immunoprecipitation IHC—Immunohistochemistry IC—Immunocytochemistry IF—Immunofluorescence
Species Cross-Reactivity Key: H—human M—mouse R—rat Hm—hamster Mk—monkey Mi—mink C—chicken X—Xenopus
Species enclosed in parentheses are predicted to react based on 100% sequence homology.

F—Flow cytometry E—ELISA D—DELFIAP®
Z—zebra fish B—bovine All—all species expected

Western Immunoblotting Protocol (Primary Antibody Incubation in BSA)

For Western blots, incubate membrane with diluted antibody in 5% w/v BSA, 1X TBS, 0.1% Tween-20 at 4°C with gentle shaking, overnight.

A Solutions and Reagents

NOTE: Prepare solutions with Milli-Q or equivalently purified water.

- 1X Phosphate Buffered Saline (PBS)
- 1X SDS Sample Buffer:** 62.5 mM Tris-HCl (pH 6.8 at 25°C), 2% w/v SDS, 10% glycerol, 50 mM DTT, 0.01% w/v bromophenol blue or phenol red
- Transfer Buffer:** 25 mM Tris base, 0.2 M glycine, 20% methanol (pH 8.5)
- 10X Tris Buffered Saline (TBS):** To prepare 1 liter of 10X TBS: 24.2 g Tris base, 80 g NaCl; adjust pH to 7.6 with HCl (use at 1X).
- Nonfat Dry Milk (weight to volume [w/v])
- Blocking Buffer:** 1X TBS, 0.1% Tween-20 with 5% w/v nonfat dry milk; for 150 ml, add 15 ml 10X TBS to 135 ml water, mix. Add 7.5 g nonfat dry milk and mix well. While stirring, add 0.15 ml Tween-20 (100%).
- Wash Buffer:** 1X TBS, 0.1% Tween-20 (TBS/T)
- Bovine Serum Albumin (BSA)
- Primary Antibody Dilution Buffer:** 1X TBS, 0.1% Tween-20 with 5% BSA; for 20 ml, add 2 ml 10X TBS to 18 ml water, mix. Add 1.0 g BSA and mix well. While stirring, add 20 µl Tween-20 (100%).
- Phototope[®]-HRP Western Blot Detection System #7071:** Includes biotinylated protein ladder, secondary anti-rabbit (#7074) antibody conjugated to horseradish peroxidase (HRP), anti-biotin antibody conjugated to HRP, LumiGLO[®] chemiluminescent reagent and peroxide.
- Prestained Protein Marker, Broad Range (Premixed Format) #7720
- Biotinylated Protein Ladder Detection Pack #7727
- Blotting Membrane:** This protocol has been optimized for nitrocellulose membranes, which CST recommends. PVDF membranes may also be used.

B Protein Blotting

A general protocol for sample preparation is described below.

- Treat cells by adding fresh media containing regulator for desired time.
- Aspirate media from cultures; wash cells with 1X PBS; aspirate.
- Lyse cells by adding 1X SDS sample buffer (100 µl per well of 6-well plate or 500 µl per plate of 10 cm diameter plate). Immediately scrape the cells off the plate and transfer the extract to a microcentrifuge tube. Keep on ice.
- Sonicate for 10–15 seconds to shear DNA and reduce sample viscosity.
- Heat a 20 µl sample to 95–100°C for 5 minutes; cool on ice.
- Microcentrifuge for 5 minutes.
- Load 20 µl onto SDS-PAGE gel (10 cm x 10 cm).

NOTE: CST recommends loading prestained molecular weight markers (#7720, 10 µl/lane) to verify electrotransfer and biotinylated protein ladder (#7727, 10 µl/lane) to determine molecular weights.

- Electrotransfer to nitrocellulose or PVDF membrane.

C Membrane Blocking and Antibody Incubations

NOTE: Volumes are for 10 cm x 10 cm (100 cm²) of membrane; for different sized membranes, adjust volumes accordingly.

- (Optional) After transfer, wash nitrocellulose membrane with 25 ml TBS for 5 minutes at room temperature.
- Incubate membrane in 25 ml of blocking buffer for 1 hour at room temperature.
- Wash three times for 5 minutes each with 15 ml of TBS/T.
- Incubate membrane and primary antibody (at the appropriate dilution) in 10 ml primary antibody dilution buffer with gentle agitation overnight at 4°C.
- Wash three times for 5 minutes each with 15 ml of TBS/T.
- Incubate membrane with HRP-conjugated secondary antibody (1:2000) and HRP-conjugated anti-biotin antibody (1:1000) to detect biotinylated protein markers in 10 ml of blocking buffer with gentle agitation for 1 hour at room temperature.
- Wash three times for 5 minutes each with 15 ml of TBS/T.

D Detection of Proteins

- Incubate membrane with 10 ml LumiGLO[®] (0.5 ml 20X LumiGLO[®], 0.5 ml 20X Peroxide and 9.0 ml Milli-Q water) with gentle agitation for 1 minute at room temperature.

NOTE: LumiGLO[®] substrate can be further diluted if signal response is too fast.

- Drain membrane of excess developing solution (do not let dry), wrap in plastic wrap and expose to x-ray film. An initial 10-second exposure should indicate the proper exposure time.

NOTE: Due to the kinetics of the detection reaction, signal is most intense immediately following LumiGLO[®] incubation and declines over the following 2 hours.

Immunofluorescence Protocol

***IMPORTANT:** Please refer to the **APPLICATIONS** section on the front page of the data sheet to determine **IF THIS PRODUCT** is validated and approved for the specific protocol you will be using.

A Solutions and Reagents

NOTE: Prepare solutions with Milli-Q or equivalently purified water.

- 1. 10X Phosphate Buffered Saline (PBS):** To prepare 1 L add 80 g sodium chloride (NaCl), 2 g potassium chloride (KCl), 14.4 g sodium phosphate, dibasic (Na_2HPO_4) and 2.4 g potassium phosphate, monobasic (KH_2PO_4) to 1 L dH_2O . Adjust pH to 7.4.
- Formaldehyde, 16%, methanol free, Polysciences, Inc. (cat# 18814), use fresh, store opened vials at 4°C in dark, dilute in PBS for use.
- Xylene
- Ethanol, anhydrous denatured, histological grade, 100% and 95%
- Distilled water (dH_2O)
- 1X PBS/0.3% Triton X-100 (PBS/Triton):** To prepare 1 L, add 100 ml 10X PBS to 900 ml dH_2O . Add 3 ml Triton X-100 and mix.
- 10 mM Sodium Citrate Buffer:** To prepare 1 L, add 2.94 g sodium citrate trisodium salt dihydrate ($\text{C}_6\text{H}_5\text{Na}_3\text{O}_7 \cdot 2\text{H}_2\text{O}$) to 1 L dH_2O . Adjust pH to 6.0.
- 1X PBS, high salt (0.4M) (high salt PBS):** To prepare 1L, add 100 ml 10X PBS to 900 ml dH_2O . Add 23.38 g NaCl and mix.
- Fluorochrome-conjugated secondary antibody

NOTE: When using any primary or fluorochrome-conjugated secondary antibody for the first time, titrate the antibody to determine which dilution allows for the strongest specific signal with the least background for your sample.

- Prolong® Gold Antifade Reagent (Invitrogen, Eugene, OR, Cat# P36930)

B Specimen Preparation

I. Cultured Cell Lines (IF-IC)

IMPORTANT: Please check the **APPLICATIONS** section of the data sheet to verify that this product is validated and approved for **(IF-IC)**.

NOTE: This general fixation protocol will work with most antibodies and cell lines. However, we recommend you try different IF/IC fixation methods (methanol or acetone alone, aldehyde alone, or combinations of these) to identify the optimal fixation protocol for each antibody and/or cell line.

NOTE: Cells should be grown, treated, fixed, and stained directly in multiwell plates, chamber slides, or on coverslips.

- Rinse cells briefly in PBS.
- Aspirate PBS, cover cells to a depth of 2-3 mm with 2-4% formaldehyde in PBS.

NOTE: Formaldehyde is toxic, use only in fume hood.

- Allow cells to fix for 15 minutes at room temperature.
- Aspirate fixative, rinse three times in PBS for 5 minutes each.

5. Methanol Permeabilization Step (if required, please refer to front page): After formaldehyde fixation, cover cells with ice-cold 100% methanol (use enough to cover cells completely to a depth of 3-5 mm, **DO NOT LET CELLS DRY**), incubate cells in methanol for 10 minutes in freezer, rinse in PBS for 5 minutes.

- Proceed with Immunostaining section C.

II. Paraffin Sections (IF-P)

IMPORTANT: Please check the **APPLICATIONS** section of the data sheet to verify that this product is validated and approved for **(IF-P)**.

Deparaffinization/Rehydration:

- Incubate sections in three washes of xylene for 5 minutes each.
- Incubate sections in two washes of 100% ethanol for 10 minutes each.
- Incubate sections in two washes of 95% ethanol for 10 minutes each.
- Rinse sections twice in dH_2O for 5 minutes each.

Antigen Unmasking:

- Place slides in room temperature 10 mM sodium citrate buffer pH 6.0.
- Bring slides to boiling in sodium citrate buffer using water bath or microwave, then maintain at 95-99°C for 10 minutes.
- Cool slides for 30 minutes on bench top.
- Rinse sections in dH_2O three times for 5 minutes each.
- Rinse sections in PBS for 5 minutes.
- Proceed with Immunostaining section C.

III. Frozen/Cryostat Sections (IF-F)

IMPORTANT: Please check the **APPLICATIONS** section of the data sheet to verify that this product is validated and approved for **(IF-F)**.

NOTE: Fresh frozen/unfixed sections should be fixed immediately in 2-4% formaldehyde as follows to preserve signaling epitopes.

- Cover sections with 2-4% formaldehyde in PBS

NOTE: Formaldehyde is toxic, use only in fume hood.

- Allow cells to fix for 15 minutes at room temperature.
- Rinse slides three times in PBS for 5 minutes each.

C Immunostaining

NOTE: All subsequent incubations should be carried out at room temperature unless otherwise noted in a humid light-tight box or covered dish/plate to prevent drying and fluorochrome fading.

- Block specimen in 5% normal serum from same species as secondary antibody (eg. normal goat serum, normal donkey serum) in PBS/Triton for 60 minutes.
- While blocking, prepare primary antibody by diluting as indicated on datasheet in PBS/Triton. You will need 50-100 μl per section, 25-50 μl per coverslip, chamber, or well (48 or 96 well plate).
- Aspirate blocking solution, apply diluted primary antibody.

NOTE: For double-labeling, prepare a cocktail of mouse and rabbit primary antibodies at their appropriate dilutions in PBS/Triton.

- Incubate **overnight** at 4°C.
- Rinse three times in PBS for 5 minutes each.

OPTION: To decrease background stain, rinse in high salt PBS for two minutes between second and third PBS rinses. Be aware, this may reduce specific staining of some antibodies.

NOTE: If using primary antibodies directly conjugated with AlexaFluor® fluorochromes, then skip to step C8.

- Incubate in fluorochrome-conjugated secondary antibody diluted in PBS/Triton for 1-2 hours at room temperature in dark.

NOTE: For double-labeling, prepare a cocktail of fluorochrome-conjugated anti-mouse and anti-rabbit primary antibodies at their appropriate dilutions in PBS/Triton.

- Rinse in PBS/high salt PBS as in step 5.
- Coverslip slides with Prolong® Gold Antifade Reagent or apply just enough to cover cells in multiwell plate.
- Seal slides by painting around edges of coverslips with nail polish.
- Examine specimens immediately using appropriate excitation wavelength, depending on fluorochrome for best results or store flat at 4°C in dark.

Immunoprecipitation Protocol / (For Analysis By Western Immunoblotting)

A Solutions and Reagents

NOTE: Prepare solutions with Milli-Q or equivalently purified water.

- 1X Phosphate Buffered Saline (PBS)
- 1X Cell Lysis Buffer:** 20 mM Tris (pH 7.5), 150 mM NaCl, 1 mM EDTA, 1 mM EGTA, 1% Triton X-100, 2.5 mM Sodium pyrophosphate, 1 mM β -glycerophosphate, 1 mM Na_3VO_4 , 1 $\mu\text{g}/\text{ml}$ Leupeptin

NOTE: Add 1 mM PMSF immediately prior to use.

- Transfer Buffer:** 25 mM Tris base, 0.2 mM glycine, 20% methanol (pH 8.5)
- Protein A or G Agarose Beads:** (Can be stored for 2 weeks at 4°C.) Please prepare according to manufacturer's instructions. Use Protein A for rabbit IgG pull down and Protein G for mouse IgG pull down.
- 3X SDS Sample Buffer:** 187.5 mM Tris-HCl (pH 6.8 at 25°C), 6% w/v SDS, 30% glycerol, 150 mM DTT, 0.03% w/v bromophenol blue

B Preparing Cell Lysates

- Aspirate media. Treat cells by adding fresh media containing regulator for desired time.
- To harvest cells under nondenaturing conditions, remove media and rinse cells once with ice-cold PBS.
- Remove PBS and add 0.5 ml ice-cold 1X cell lysis buffer to each plate (10 cm) and incubate the plates on ice for 5 minutes.
- Scrape cells off the plates and transfer to microcentrifuge tubes. Keep on ice.
- Sonicate samples on ice three times for 5 seconds each.
- Microcentrifuge for 10 minutes at 14,000 X g, 4°C, and transfer the supernatant to a new tube. If necessary, lysate can be stored at -80°C.

C Immunoprecipitation

Optional: It may be necessary to perform a lysate pre-clearing step to reduce non-specific binding to the Protein A/G agarose beads (See section below).

- Take 200 μl cell lysate and add primary antibody. Incubate with gentle rocking overnight at 4°C.
- Add either protein A or G agarose beads (20 μl of 50% bead slurry). Incubate with gentle rocking for 1–3 hours at 4°C.
- Microcentrifuge for 30 seconds at 4°C. Wash pellet five times with 500 μl of 1X cell lysis buffer. Keep on ice during washes.
- Resuspend the pellet with 20 μl 3X SDS sample buffer. Vortex, then microcentrifuge for 30 seconds.
- Heat the sample to 95–100°C for 2–5 minutes and microcentrifuge for 1 minute at 14,000 X g.
- Load the sample (15–30 μl) on SDS-PAGE gel (12–15%).
- Analyze sample by Western blotting (see Western Immunoblotting Protocol).

Cell Lysate Pre-Clearing (Optional)

- Take 200 μl cell lysate and add to either Protein A or G agarose beads (20 μl of 50% bead slurry).
- Incubate at 4°C for 30 – 60 minutes.
- Spin for 10 minutes at 4°C. Transfer the supernatant to a fresh tube.
- Proceed to step 1 of Immunoprecipitation.

Flow Cytometry Protocol for Intracellular Staining Using Conjugated Secondary Antibodies

A Solutions and Reagents

- 1X Phosphate Buffered Saline (PBS):** Dissolve 8 g NaCl, 0.2 g KCl, 1.44 g Na_2HPO_4 and 0.24 g KH_2PO_4 in 800 mL distilled water (dH_2O). Adjust the pH to 7.4 with HCl and the volume to 1 liter. Store at room temperature.
- Formaldehyde (methanol free)
- Incubation Buffer:** Dissolve 0.5 g bovine serum albumin (BSA) in 100mL 1X PBS. Store at 4°C

B Fixation

- Collect cells by centrifugation and aspirate supernatant.
- Resuspend cells briefly in 0.5-1 ml PBS. Add formaldehyde to a final concentration of 2-4% formaldehyde.
- Fix for 10 minutes at 37°C.
- Chill tubes on ice for 1 minute.

C Permeabilization

- Permeabilize cells by adding ice-cold 100% methanol slowly to pre-chilled cells, while gently vortexing, to a final concentration of 90% methanol. Alternatively, to remove fix prior to permeabilization, pellet cells by centrifugation and resuspend in 90% methanol.
- Incubate 30 minutes on ice.
- Proceed with staining or store cells at -20°C in 90% methanol.

D Staining Using Unlabeled Primary and Conjugated Secondary Antibodies

NOTE: Allow for isotype matched controls for monoclonal antibodies or species matched IgG for polyclonal antibodies. Count cells using a hemacytometer or alternative method.

- Aliquot 0.5-1x10⁶ cells into each assay tube (by volume).
- Add 2-3 ml Incubation Buffer to each tube and rinse by centrifugation. Repeat.
- Resuspend cells in 100 μl Incubation Buffer per assay tube.
- Block in Incubation Buffer for 10 minutes at room temperature.
- Add the primary antibody at the appropriate dilution to the assay tubes (see individual antibody data sheet for the appropriate dilution).
- Incubate for 30-60 minutes at room temperature.
- Rinse as before in Incubation Buffer by centrifugation.
- Resuspend cells in fluorochrome-conjugated secondary antibody*, diluted in Incubation Buffer according to the manufacturer's recommendations.
- Incubate for 30 minutes at room temperature.
- Rinse as before in Incubation Buffer by centrifugation.
- Resuspend cells in 0.5 ml PBS and analyze on flow cytometer.

*Recommended Secondary Antibodies from Invitrogen.

A-11070 Alexa Fluor® 488 F(ab')₂ fragment of goat anti-rabbit IgG (H+L) (1:1000 dilution)
A-11017 Alexa Fluor® 488 F(ab')₂ fragment of goat anti-mouse IgG (H+L) (1:1000 dilution)