

ASF1B (C70E2) Rabbit mAb

100 µl
 (10 Western mini-blots)

Orders ■ 877-616-CELL (2355)
 orders@cellsignal.com
Support ■ 877-678-TECH (8324)
 info@cellsignal.com
Web ■ www.cellsignal.com

New 05/08


Products are sold and licensed for use by end-users only for research or laboratory use. Products are not to be used on humans for clinical diagnostic, drug screening or testing purposes.

Applications	Species Cross-Reactivity*	Molecular Wt.	Source	Isotype
W, IP, IF-IC Endogenous	H, Mk	19 kDa	Rabbit**	IgG


Background: ASF1 was first identified in *S. cerevisiae* based on its ability to de-repress transcriptional silencing when overexpressed (1). While only one gene exists in yeast and *Drosophila*, mammalian cells contain the two highly homologous ASF1A and ASF1B genes (2). ASF1A and ASF1B function as histone chaperones, delivering histone H3/H4 dimers to CAF-1 or HIRA histone deposition complexes to facilitate replication-coupled and replication-independent nucleosome assembly on DNA (2-5). Both ASF1A and ASF1B bind to CAF-1, but only ASF1A binds to HIRA (5). In addition to playing a role in DNA replication and gene silencing, ASF1 functions in DNA damage repair, genome stability and cellular senescence. Deletion of ASF1 in yeast and *Drosophila* confers sensitivity to various DNA damaging agents and inhibitors of DNA replication, increases genomic instability and sister chromatid exchange, and activates the DNA damage checkpoint (6-8). Depletion of both ASF1A and ASF1B in mammalian cells results in the accumulation of cells in S phase, increased phosphorylation of H2A.X, centrosome amplification and apoptosis (9,10). ASF1A is required for the formation of senescence-associated heterochromatin foci (SAHF), with overexpression of ASF1A inducing senescence in primary cells (4). Both ASF1A and ASF1B are phosphorylated in S phase by the Tousled-like kinases TLK1 and TLK2, and are dephosphorylated when TLK1 and TLK2 are inactivated by Chk1 kinase in response to replicative stress (11,12). The function of ASF1 phosphorylation is not yet understood.

Specificity/Sensitivity: ASF1B (C70E2) Rabbit mAb detects endogenous levels of total ASF1B protein. The antibody does not cross-react with ASF1A protein.

Source/Purification: Monoclonal antibody is produced by immunizing animals with a synthetic peptide (KLH-coupled) corresponding to the carboxy-terminus of the human ASF1B protein.

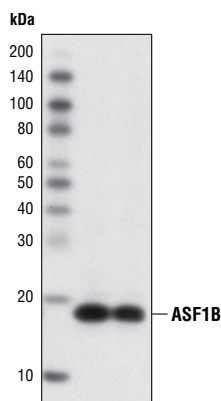

Western blot analysis of extracts from HeLa and HCT116 cells using ASF1B (C70E2) Rabbit mAb.

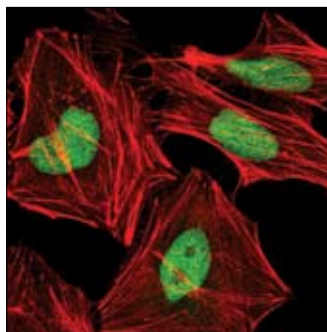

Confocal immunofluorescent analysis of HeLa cells using ASF1B (C70E2) Rabbit mAb (green). Actin filaments have been labeled with DY-554 phalloidin (red).

Entrez-Gene ID #55723
Swiss-Prot Acc. #Q53G51

Storage: Supplied in 10 mM sodium HEPES (pH 7.5), 150 mM NaCl, 100 µg/ml BSA, 50% glycerol and less than 0.02% sodium azide. Store at -20°C. Do not aliquot the antibody.

***Species cross-reactivity is determined by Western blot.**

****Anti-rabbit secondary antibodies must be used to detect this antibody.**

Recommended Antibody Dilutions:

Western blotting	1:1000
Immunoprecipitation	1:100
Immunofluorescence (IF-IC)	1:100

Companion Products:

- Phospho-TLK1 (Ser695) Antibody #4121
- TLK1 Antibody #4125
- Phospho-Histone H2A.X (Ser139) Antibody #2577
- Phospho-Histone H2A.X (Ser139) (20E3) Rabbit mAb #9718
- Histone H2A.X Antibody #2595
- Phototope®-HRP Western Blot Detection System, Anti-rabbit IgG, HRP-linked Antibody #7071
- Anti-rabbit IgG, HRP-linked Antibody #7074
- Prestained Protein Marker, Broad Range (Premixed Format) #7720
- Biotinylated Protein Ladder Detection Pack #7727
- 20X LumiGLO® Reagent and 20X Peroxide #7003

Background References:

- (1) Singer, M.S. et al. (1998) *Genetics* 150, 613-632.
- (2) Mousson, F. et al. (2006) *Chromosoma* 116, 79-93.
- (3) Tang, Y. et al. (2006) *Nat. Struct. Mol. Biol.* 13, 921-929.
- (4) Zhang, R. et al. (2005) *Dev. Cell.* 8, 19-30.
- (5) Daganzo, S.M. et al. (2003) *Curr. Biol.* 13, 2148-2158.
- (6) Ramey, C.J. et al. (2004) *Mol. Cell. Biol.* 24, 10313-10327.
- (7) Prado, F. et al. (2004) *EMBO Rep.* 5, 497-502.
- (8) Tyler, J.K. et al. (1999) *Nature* 402, 555-560.
- (9) Sanematsu, F. et al. (2006) *J. Biol. Chem.* 281, 13817-13827.
- (10) Groth, A. et al. (2005) *Mol. Cell.* 17, 301-311.
- (11) Silljé, H.H. and Nigg, E.A. (2001) *Curr. Biol.* 11, 1068-1073.
- (12) Carrera, P. et al. (2003) *Genes Dev.* 17, 2578-2590.

IMPORTANT: For Western blots, incubate membrane with diluted antibody in 5% w/v BSA, 1X TBS, 0.1% Tween-20 at 4°C with gentle shaking, overnight.

Western Immunoblotting Protocol (Primary Antibody Incubation in BSA)

For Western blots, incubate membrane with diluted antibody in 5% w/v BSA, 1X TBS, 0.1% Tween-20 at 4°C with gentle shaking, overnight.

A Solutions and Reagents

NOTE: Prepare solutions with Milli-Q or equivalently purified water.

- 1X Phosphate Buffered Saline (PBS)
- 1X SDS Sample Buffer:** 62.5 mM Tris-HCl (pH 6.8 at 25°C), 2% w/v SDS, 10% glycerol, 50 mM DTT, 0.01% w/v bromophenol blue or phenol red
- Transfer Buffer:** 25 mM Tris base, 0.2 M glycine, 20% methanol (pH 8.5)
- 10X Tris Buffered Saline (TBS):** To prepare 1 liter of 10X TBS: 24.2 g Tris base, 80 g NaCl; adjust pH to 7.6 with HCl (use at 1X).
- Nonfat Dry Milk (weight to volume [w/v])
- Blocking Buffer:** 1X TBS, 0.1% Tween-20 with 5% w/v nonfat dry milk; for 150 ml, add 15 ml 10X TBS to 135 ml water, mix. Add 7.5 g nonfat dry milk and mix well. While stirring, add 0.15 ml Tween-20 (100%).
- Wash Buffer:** 1X TBS, 0.1% Tween-20 (TBS/T)
- Bovine Serum Albumin (BSA)
- Primary Antibody Dilution Buffer:** 1X TBS, 0.1% Tween-20 with 5% BSA; for 20 ml, add 2 ml 10X TBS to 18 ml water, mix. Add 1.0 g BSA and mix well. While stirring, add 20 µl Tween-20 (100%).
- Phototope®-HRP Western Blot Detection System #7071:** Includes biotinylated protein ladder, secondary anti-rabbit (#7074) antibody conjugated to horseradish peroxidase (HRP), anti-biotin antibody conjugated to HRP, LumiGLO® chemiluminescent reagent and peroxide.
- Prestained Protein Marker, Broad Range (Premixed Format) #7720
- Biotinylated Protein Ladder Detection Pack #7727
- Blotting Membrane:** This protocol has been optimized for nitrocellulose membranes, which CST recommends. PVDF membranes may also be used.

B Protein Blotting

A general protocol for sample preparation is described below.

- Treat cells by adding fresh media containing regulator for desired time.
- Aspirate media from cultures; wash cells with 1X PBS; aspirate.
- Lyse cells by adding 1X SDS sample buffer (100 µl per well of 6-well plate or 500 µl per plate of 10 cm diameter plate). Immediately scrape the cells off the plate and transfer the extract to a microcentrifuge tube. Keep on ice.
- Sonicate for 10–15 seconds to shear DNA and reduce sample viscosity.
- Heat a 20 µl sample to 95–100°C for 5 minutes; cool on ice.
- Microcentrifuge for 5 minutes.
- Load 20 µl onto SDS-PAGE gel (10 cm x 10 cm).

NOTE: CST recommends loading prestained molecular weight markers (#7720, 10 µl/lane) to verify electrotransfer and biotinylated protein ladder (#7727, 10 µl/lane) to determine molecular weights.

- Electrotransfer to nitrocellulose or PVDF membrane.

C Membrane Blocking and Antibody Incubations

NOTE: Volumes are for 10 cm x 10 cm (100 cm²) of membrane; for different sized membranes, adjust volumes accordingly.

- (Optional) After transfer, wash nitrocellulose membrane with 25 ml TBS for 5 minutes at room temperature.
- Incubate membrane in 25 ml of blocking buffer for 1 hour at room temperature.
- Wash three times for 5 minutes each with 15 ml of TBS/T.
- Incubate membrane and primary antibody (at the appropriate dilution) in 10 ml primary antibody dilution buffer with gentle agitation overnight at 4°C.
- Wash three times for 5 minutes each with 15 ml of TBS/T.
- Incubate membrane with HRP-conjugated secondary antibody (1:2000) and HRP-conjugated anti-biotin antibody (1:1000) to detect biotinylated protein markers in 10 ml of blocking buffer with gentle agitation for 1 hour at room temperature.
- Wash three times for 5 minutes each with 15 ml of TBS/T.

D Detection of Proteins

- Incubate membrane with 10 ml LumiGLO® (0.5 ml 20X LumiGLO®, 0.5 ml 20X Peroxide and 9.0 ml Milli-Q water) with gentle agitation for 1 minute at room temperature.

NOTE: LumiGLO® substrate can be further diluted if signal response is too fast.

- Drain membrane of excess developing solution (do not let dry), wrap in plastic wrap and expose to x-ray film. An initial 10-second exposure should indicate the proper exposure time.

NOTE: Due to the kinetics of the detection reaction, signal is most intense immediately following LumiGLO® incubation and declines over the following 2 hours.

Immunoprecipitation Protocol / (For Analysis By Western Immunoblotting)

A Solutions and Reagents

NOTE: Prepare solutions with Milli-Q or equivalently purified water.

1. 1X Phosphate Buffered Saline (PBS)
2. **1X Cell Lysis Buffer:** 20 mM Tris (pH 7.5), 150 mM NaCl, 1 mM EDTA, 1 mM EGTA, 1% Triton X-100, 2.5 mM Sodium pyrophosphate, 1 mM β -glycerophosphate, 1 mM Na_3VO_4 , 1 $\mu\text{g/ml}$ Leupeptin

NOTE: Add 1 mM PMSF immediately prior to use.

1. **Transfer Buffer:** 25 mM Tris base, 0.2 mM glycine, 20% methanol (pH 8.5)
2. **Protein A or G Agarose Beads:** (Can be stored for 2 weeks at 4°C.) Please prepare according to manufacturer's instructions. Use Protein A for rabbit IgG pull down and Protein G for mouse IgG pull down.
3. **3X SDS Sample Buffer:** 187.5 mM Tris-HCl (pH 6.8 at 25°C), 6% w/v SDS, 30% glycerol, 150 mM DTT, 0.03% w/v bromophenol blue

B Preparing Cell Lysates

1. Aspirate media. Treat cells by adding fresh media containing regulator for desired time.
2. To harvest cells under nondenaturing conditions, remove media and rinse cells once with ice-cold PBS.
3. Remove PBS and add 0.5 ml ice-cold 1X cell lysis buffer to each plate (10 cm) and incubate the plates on ice for 5 minutes.
4. Scrape cells off the plates and transfer to microcentrifuge tubes. Keep on ice.
5. Sonicate samples on ice three times for 5 seconds each.
6. Microcentrifuge for 10 minutes at 14,000 X g, 4°C, and transfer the supernatant to a new tube. If necessary, lysate can be stored at -80°C.

C Immunoprecipitation

Optional: It may be necessary to perform a lysate pre-clearing step to reduce non-specific binding to the Protein A/G agarose beads (See section below).

1. Take 200 μl cell lysate and add primary antibody. Incubate with gentle rocking overnight at 4°C.
2. Add either protein A or G agarose beads (20 μl of 50% bead slurry). Incubate with gentle rocking for 1–3 hours at 4°C.
3. Microcentrifuge for 30 seconds at 4°C. Wash pellet five times with 500 μl of 1X cell lysis buffer. Keep on ice during washes.
4. Resuspend the pellet with 20 μl 3X SDS sample buffer. Vortex, then microcentrifuge for 30 seconds.
5. Heat the sample to 95–100°C for 2–5 minutes and microcentrifuge for 1 minute at 14,000 X g.
6. Load the sample (15–30 μl) on SDS-PAGE gel (12–15%).
7. Analyze sample by Western blotting (see Western Immunoblotting Protocol).

Cell Lysate Pre-Clearing (Optional)

1. Take 200 μl cell lysate and add to either Protein A or G agarose beads (20 μl of 50% bead slurry).
2. Incubate at 4°C for 30 – 60 minutes.
3. Spin for 10 minutes at 4°C. Transfer the supernatant to a fresh tube.
4. Proceed to step 1 of Immunoprecipitation.

Immunofluorescence Protocol

***IMPORTANT:** Please refer to the **APPLICATIONS** section on the front page of the data sheet to determine **IF THIS PRODUCT** is validated and approved for the specific protocol you will be using.

A Solutions and Reagents

NOTE: Prepare solutions with Milli-Q or equivalently purified water.

- 1. 10X Phosphate Buffered Saline (PBS):** To prepare 1 L add 80 g sodium chloride (NaCl), 2 g potassium chloride (KCl), 14.4 g sodium phosphate, dibasic (Na_2HPO_4) and 2.4 g potassium phosphate, monobasic (KH_2PO_4) to 1 L dH_2O . Adjust pH to 7.4.
- Formaldehyde, 16%, methanol free, Polysciences, Inc. (cat# 18814), use fresh, store opened vials at 4°C in dark, dilute in PBS for use.
- Xylene
- Ethanol, anhydrous denatured, histological grade, 100% and 95%
- Distilled water (dH_2O)
- Blocking Buffer:** To prepare 25 mL, add 2.5 mL 10X PBS, 1.25 mL normal serum from the same species as the secondary antibody (eg. normal goat serum, normal donkey serum) and 21.25 mL dH_2O and mix well. While stirring, add 75 μL Triton X-100 (100%).
- Antibody Dilution Buffer:** To prepare 40 mL, add 4 mL 10X PBS to 36 mL dH_2O , mix. Add 0.4 g BSA and mix well. While stirring, add 120 μL Triton X-100 (100%).
- 10 mM Sodium Citrate Buffer:** To prepare 1 L, add 2.94 g sodium citrate trisodium salt dihydrate ($\text{C}_6\text{H}_5\text{Na}_3\text{O}_7 \cdot 2\text{H}_2\text{O}$) to 1 L dH_2O . Adjust pH to 6.0.
- 1X PBS, high salt (0.4M) (high salt PBS):** To prepare 1L, add 100 ml 10X PBS to 900 ml dH_2O . Add 23.38 g NaCl and mix.
- Fluorochrome-conjugated secondary antibody

NOTE: When using any primary or fluorochrome-conjugated secondary antibody for the first time, titrate the antibody to determine which dilution allows for the strongest specific signal with the least background for your sample.

- Prolong® Gold Antifade Reagent (Invitrogen, Eugene, OR, Cat# P36930)

B Specimen Preparation

I. Cultured Cell Lines (IF-IC)

IMPORTANT: Please check the **APPLICATIONS** section of the data sheet to verify that this product is validated and approved for **(IF-IC)**.

NOTE: Cells should be grown, treated, fixed, and stained directly in multiwell plates, chamber slides, or on coverslips.

- Rinse cells briefly in PBS.
- Aspirate PBS, cover cells to a depth of 2-3 mm with 2-4% formaldehyde in PBS.

NOTE: Formaldehyde is toxic, use only in fume hood.

- Allow cells to fix for 15 minutes at room temperature.
- Aspirate fixative, rinse three times in PBS for 5 minutes each.
- Methanol Permeabilization Step (if required, please refer to front page):** After formaldehyde fixation, cover cells with ice-cold 100% methanol (use enough to cover cells completely to a depth of 3-5 mm, DO NOT LET CELLS DRY), incubate cells in methanol for 10 minutes at -20°C , rinse in PBS for 5 minutes.
- Proceed with Immunostaining section C.

II. Paraffin Sections (IF-P)

IMPORTANT: Please check the **APPLICATIONS** section of the data sheet to verify that this product is validated and approved for **(IF-P)**.

Deparaffinization/Rehydration:

- Incubate sections in three washes of xylene for 5 minutes each.
- Incubate sections in two washes of 100% ethanol for 10 minutes each.
- Incubate sections in two washes of 95% ethanol for 10 minutes each.
- Rinse sections twice in dH_2O for 5 minutes each.

Antigen Unmasking:

- Place slides in room temperature 10 mM sodium citrate buffer pH 6.0.
- Bring slides to boiling in sodium citrate buffer using water bath or microwave, then maintain at $95-99^\circ\text{C}$ for 10 minutes.
- Cool slides for 30 minutes on bench top.
- Rinse sections in dH_2O three times for 5 minutes each.
- Rinse sections in PBS for 5 minutes.
- Proceed with Immunostaining section C.

III. Frozen/Cryostat Sections (IF-F)

IMPORTANT: Please check the **APPLICATIONS** section of the data sheet to verify that this product is validated and approved for **(IF-F)**.

NOTE: Fresh frozen/unfixed sections should be fixed immediately in 2-4% formaldehyde as follows to preserve signaling epitopes.

- Cover sections with 2-4% formaldehyde in PBS

NOTE: Formaldehyde is toxic, use only in fume hood.

- Allow sections to fix for 15 minutes at room temperature.
- Rinse slides three times in PBS for 5 minutes each.

C Immunostaining

NOTE: All subsequent incubations should be carried out at room temperature unless otherwise noted in a humid light-tight box or covered dish/plate to prevent drying and fluorochrome fading.

- Block specimen in Blocking Buffer for 60 minutes.
- While blocking, prepare primary antibody by diluting as indicated on datasheet in Antibody Dilution Buffer.
- Aspirate blocking solution, apply diluted primary antibody.

NOTE: For double-labeling, prepare a cocktail of the primary antibodies at their appropriate dilution in Antibody Dilution Buffer.

- Incubate overnight at 4°C .
- Rinse three times in PBS for 5 minutes each.

OPTION: To decrease background stain, rinse in high salt PBS for two minutes between second and third PBS rinses. Be aware, this may reduce specific staining of some antibodies.

NOTE: If using primary antibodies directly conjugated with Alexa Fluor® fluorochromes, then skip to step C8.

- Incubate specimen in fluorochrome-conjugated secondary antibody diluted in Antibody Dilution Buffer for 1-2 hours at room temperature in dark.

NOTE: For double-labeling, prepare a cocktail of fluorochrome-conjugated secondary antibodies at their appropriate dilutions in Antibody Dilution Buffer.

- Rinse in PBS/high salt PBS as in step 5.
- Coverslip slides with Prolong® Gold Antifade Reagent or apply just enough to cover cells in multiwell plate.
- Seal slides by painting around edges of coverslips with nail polish.
- For best results examine specimens immediately using appropriate excitation wavelength. For long term storage, store slides flat at 4°C protected from light.